

Certificate of Analysis

Verra Lab
Verra Lab
Netherlands
support@verralab.com

Liquilabs s.r.o.
Inovační 122
252 41 Zlatníky-Hodkovice
Czechia
www.liquilabs.cz



[Verify Results Online](#)

Sample Identification

Sample Name Glow 70 mg
Batch Number
Date Published 2026-08-07 10:33

Results for LYO-0468

Peptides	Result	Unit	Uncertainty	Acceptable Range
BPC-157 Assay Peptide Screening 0.1% TFA	13.36	mg	[± 0.07]	
GHK-Cu Assay Peptide Screening 0.1% TFA	52.9	mg	[± 0.3]	
BPC-157 Purity Peptide Screening 0.1% TFA	98.8	%	[± 0.5]	
GHK-Cu Purity Peptide Screening 0.1% TFA	> 99.8	%		
BPC-157 Identification by Spectrum (FTIR) Peptide Screening 0.1% TFA	995		[± 5]	
GHK-Cu Identification by Spectrum (FTIR) Peptide Screening 0.1% TFA	989		[± 5]	
BPC-157 Identification by RT Peptide Screening 0.1% TFA	0.996		[± 0.005]	
GHK-Cu Identification by RT Peptide Screening 0.1% TFA	0.997		[± 0.005]	
Thymosin Beta 4 (TB-500) Assay Peptide Screening 0.1% TFA	11.70	mg	[± 0.06]	
Thymosin Beta 4 (TB-500) Purity Peptide Screening 0.1% TFA	99.5	%	[± 0.5]	
Thymosin Beta 4 (TB-500) Identification by Spectrum (FTIR) Peptide Screening 0.1% TFA	988		[± 5]	
Thymosin Beta 4 (TB-500) Identification by RT Peptide Screening 0.1% TFA	0.997		[± 0.005]	
Microbiology	Result	Unit	Uncertainty	Acceptable Range
Bacterial Endotoxin Chromgenic USP<85>/ Eur. Ph. 2.6.14. Bacterial Endotoxin Chromgenic Test	0.003	EU/mg	[± 0]	
Mass Spectrometry	Result	Unit	Uncertainty	Acceptable Range
Molecular Ion Mass Identification (BPC-157) Mass Spectrometry Identity	1418	Da	[± 1]	
Molecular Ion Mass Identification (GHK-Cu) Mass Spectrometry Identity	341	Da	[± 1]	
Molecular Ion Mass Identification (TB-500) Mass Spectrometry Identity	4963	Da	[± 1]	

	Method Specification	
Determination of identity, content and purity of GHK-Cu		
<i>Document number</i> GHKCu_006_2026	<i>Superseded document</i> -	<i>Number of pages</i> 4

1. Content Assessment

1.1. Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Colum Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

1.2. Chromatographic conditions

Chromatographic conditions	
Eluent A	0.1% TFA in Water (HPLC, Gradient Grade)
Eluent B	0.1% TFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.4 mL/min
Program	Gradient elution
Injection volume	2 µL
Colum Temperature	60°C
Column	Phenomenex Biozen Peptide Polar C18, 150x2.1mm 3µm
Detection wavelength	225nm


Gradient Program		
Time [min]	A [%]	B [%]
3	100	0
18	45	55
20	1	99
21	1	99
21.01	100	0
26	end	

1.3. Sample preparation

Whole amount of container was dissolved in 2mL of water (LCMS Grade). 100 µL of sample was transferred to HPLC vial and diluted by 900 µL water (LCMS Grade) and submitted for analysis.

1.4. Calibration curve

Calibration curve detail	
Quantitative method	External Standard
Calibration Type	Linear
Number of calibration points	5
Force through Zero	Enabled
Weighting Method	None

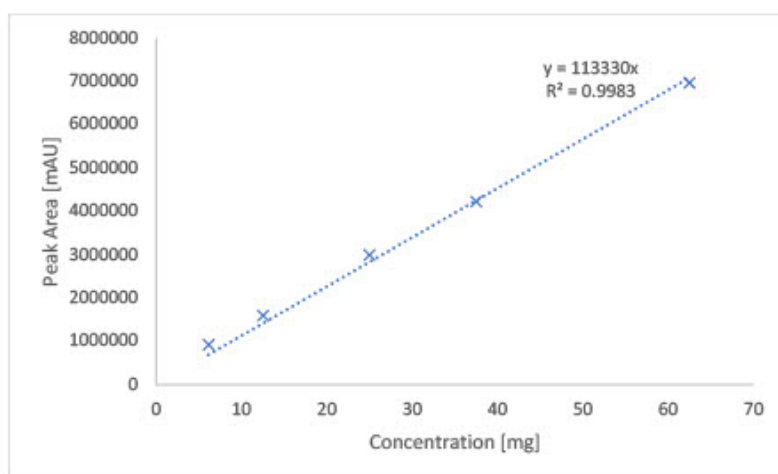

2. Purity assessment

2.1 Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Colum Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

2.2 Chromatographic conditions

Chromatographic conditions	
Eluent A	0.1% TFA in Water (HPLC, Gradient Grade)
Eluent B	0.1% TFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.4 mL/min
Program	Gradient elution
Injection volume	2 µL
Colum Temperature	60°C
Column	Phenomenex Biozen Peptide Polar C18, 150x2.1mm 3µm
Detection wavelength	225nm

Gradient Program		
Time [min]	A [%]	B [%]
3	100	0
18	45	55
20	1	99
21	1	99
21.01	100	0
26	end	

2.3 Purity assesment

Purity of compound assesed by area normalization method, comparing area of each peak to sum of area of all peaks detected at wavelenght of 225nm.

3. Identity Assessment

3.1 Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Column Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

3.2 Chromatographic conditions

Chromatographic conditions	
Eluent A	0.1% TFA in Water (HPLC, Gradient Grade)
Eluent B	0.1% TFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.4 mL/min
Program	Gradient elution
Injection volume	2 µL
Column Temperature	60°C
Column	Phenomenex Biozen Peptide Polar C18, 150x2.1mm 3µm
Mass spectrometry	Scan: positive 280-2000 Da

Gradient Program		
Time [min]	A [%]	B [%]
3	100	0
18	45	55
20	1	99
21	1	99
21.01	100	0
26	end	

3.3 Molecular Ion Mass evaluation

Molecular ion mass was determined by deconvolution of multiply charged ESI-MS spectra to calculate the average neutral (zero-charge) molecular mass by equation:

$$M(\text{neutral}) = (z_i((mz_i) - H)) - ME$$

Where:

mz_i - Measured mass of charged particle

z_i - charge

H - proton mass (1.0076 Da)

ME - mass error

	Method Specification	
Determination of identity, content and purity of BPC-157		
<i>Document number</i> BPCPol_006_2026	<i>Superseded document</i> -	<i>Number of pages</i> 4

1. Content Assessment

1.1. Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Colum Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

1.2. Chromatographic conditions

Chromatographic conditions	
Eluent A	0.1% TFA in Water (HPLC, Gradient Grade)
Eluent B	0.1% TFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.4 mL/min
Program	Gradient elution
Injection volume	2 µL
Colum Temperature	60°C
Column	Phenomenex Biozen Peptide Polar C18, 150x2.1mm 3µm
Detection wavelength	225nm


Gradient Program		
Time [min]	A [%]	B [%]
3	100	0
18	45	55
20	1	99
21	1	99
21.01	100	0
26	end	

1.3. Sample preparation

Whole amount of container was dissolved in 2mL of water (LCMS Grade). 100 µL of sample was transferred to HPLC vial and diluted by 900 µL water (LCMS Grade) and submitted for analysis.

1.4. Calibration curve

Calibration curve detail	
Quantitative method	External Standard
Calibration Type	Linear
Number of calibration points	5
Force through Zero	Enabled
Weighting Method	None

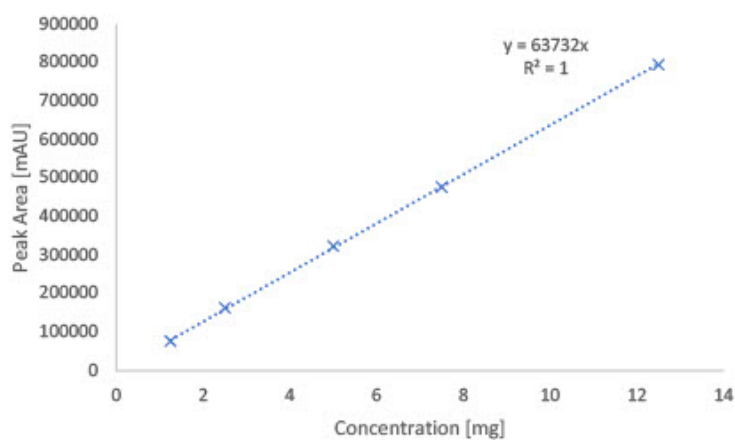

2. Purity assessment

2.1 Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Colum Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

2.2 Chromatographic conditions

Chromatographic conditions	
Eluent A	0.05% TFA in Water (HPLC, Gradient Grade)
Eluent B	0.0425% TFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.9 mL/min
Program	Gradient elution
Injection volume	2 µL
Colum Temperature	55°C
Column	Phenomenex Biozen Peptide Polar C18, 150x2.1mm 3µm
Detection wavelength	225nm

Gradient Program		
Time [min]	A [%]	B [%]
3	100	0
18	45	55
20	1	99
21	1	99
21.01	100	0
26	end	

2.3 Purity assesment

Purity of compound assesed by area normalization method, comparing area of each peak to sum of area of all peaks detected at wavelenght of 225nm.

3. Identity Assessment

3.1 Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Column Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

3.2 Chromatographic conditions

Chromatographic conditions	
Eluent A	0.05% TFA in Water (HPLC, Gradient Grade)
Eluent B	0.0425% TFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.9 mL/min
Program	Gradient elution
Injection volume	2 µL
Column Temperature	55°C
Column	Phenomenex Biozen Peptide Polar C18, 150x2.1mm 3µm
Mass spectrometry	Scan: positive 280-2000 Da

Gradient Program		
Time [min]	A [%]	B [%]
3	100	0
18	45	55
20	1	99
21	1	99
21.01	100	0
26	end	

3.3 Molecular Ion Mass evaluation

Molecular ion mass was determined by deconvolution of multiply charged ESI-MS spectra to calculate the average neutral (zero-charge) molecular mass by equation:

$$M(\text{neutral}) = (z_i((mz_i) - H)) - ME$$

Where:

mz_i - Measured mass of charged particle

z_i - charge

H - proton mass (1.0076 Da)

ME – mass error

	Method Specification	
Determination of identity, content and purity of TB-500		
<i>Document number</i> TBPol_006_2026	<i>Superseded document</i> -	<i>Number of pages</i> 4

1. Content Assessment

1.1. Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Colum Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

1.2. Chromatographic conditions

Chromatographic conditions	
Eluent A	0.1% TFA in Water (HPLC, Gradient Grade)
Eluent B	0.1% TFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.4 mL/min
Program	Gradient elution
Injection volume	2 µL
Colum Temperature	60°C
Column	Phenomenex Biozen Peptide Polar C18, 150x2.1mm 3µm
Detection wavelength	225nm


Gradient Program		
Time [min]	A [%]	B [%]
3	100	0
18	45	55
20	1	99
21	1	99
21.01	100	0
26	end	

1.3. Sample preparation

Whole amount of container was dissolved in 2mL of water (LCMS Grade). 100 µL of sample was transferred to HPLC vial and diluted by 900 µL water (LCMS Grade) and submitted for analysis.

1.4. Calibration curve

Calibration curve detail	
Quantitative method	External Standard
Calibration Type	Linear
Number of calibration points	5
Force through Zero	Enabled
Weighting Method	None

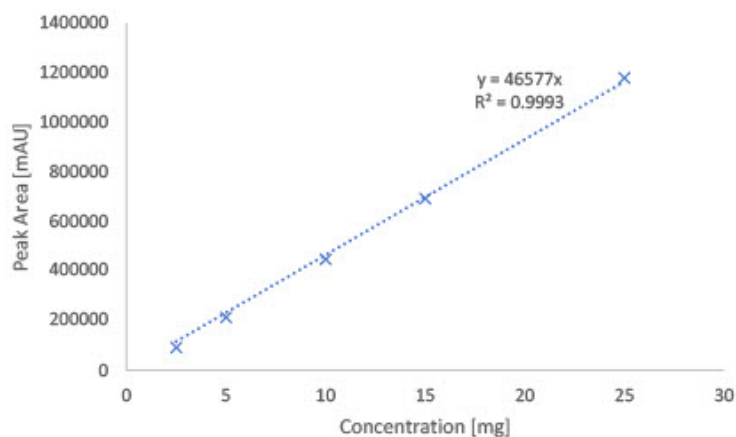

2. Purity assessment

2.1 Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Colum Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

2.2 Chromatographic conditions

Chromatographic conditions	
Eluent A	0.05% TFA in Water (HPLC, Gradient Grade)
Eluent B	0.0425% TFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.9 mL/min
Program	Gradient elution
Injection volume	2 µL
Colum Temperature	55°C
Column	Phenomenex Biozen Peptide Polar C18, 150x2.1mm 3µm
Detection wavelength	225nm

Gradient Program		
Time [min]	A [%]	B [%]
3	100	0
18	45	55
20	1	99
21	1	99
21.01	100	0
26	end	

2.3 Purity assesment

Purity of compound assesed by area normalization method, comparing area of each peak to sum of area of all peaks detected at wavelenght of 225nm.

3. Identity Assessment

3.1 Instrumentation

Module	Name	Serial Number
System Controller	Shimadzu CBM-40 Lite	L221226351398
Degassing Unit	Shimadzu DGU-403	NA
Pump	Shimadzu LC-40B XR	L22146350580
Autosampler	Shimadzu SIL-40C XR	L22216351622
Colum Thermostat	Shimadzu CTO-40S	L22236351602
PDA Detector	Shimadzu SPD-M40	L22276352808
SQ MS Detector	Shimadzu LCMS-2050	O12476200760

3.2 Chromatographic conditions

Chromatographic conditions	
Eluent A	0.05% TFA in Water (HPLC, Gradient Grade)
Eluent B	0.0425% TFA in Acetonitrile (HPLC, Gradient Grade)
Flow rate	0.9 mL/min
Program	Gradient elution
Injection volume	2 µL
Colum Temperature	55°C
Column	Phenomenex Biozen Peptide Polar C18, 150x2.1mm 3µm
Mass spectrometry	Scan: positive 280-2000 Da

Gradient Program		
Time [min]	A [%]	B [%]
3	100	0
18	45	55
20	1	99
21	1	99
21.01	100	0
26	end	

3.3 Molecular Ion Mass evaluation

Molecular ion mass was determined by deconvolution of multiply charged ESI-MS spectra to calculate the average neutral (zero-charge) molecular mass by equation:

$$M(\text{neutral}) = (z_i((mz_i) - H)) - ME$$

Where:

mz_i - Measured mass of charged particle

z_i - charge

H - proton mass (1.0076 Da)

ME – mass error

Analysis Report



LIQUILABS

Analytical services


Analysis of quantity and purity of active ingredient by UHPLC with UV detection

Sample Information


Injection Volume : 2
Data File : LYO-0468__009.lcd
Method File : Peptide screening polar_KLOW.lcm
Date Acquired : 8/5/2026 1:22:21 PM

Chromatogram

uAU

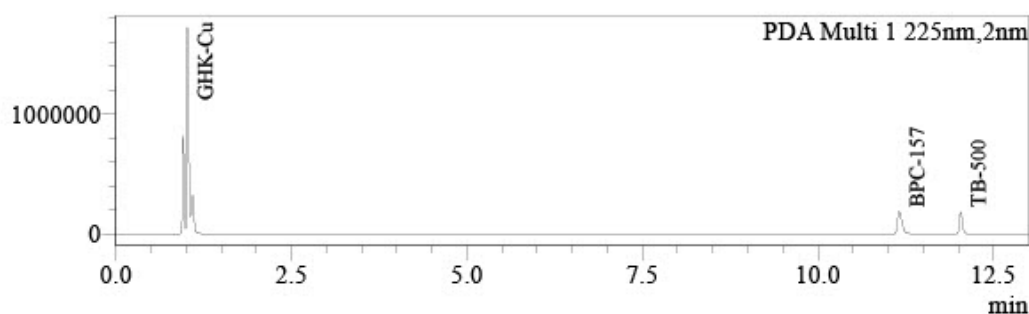

uAU

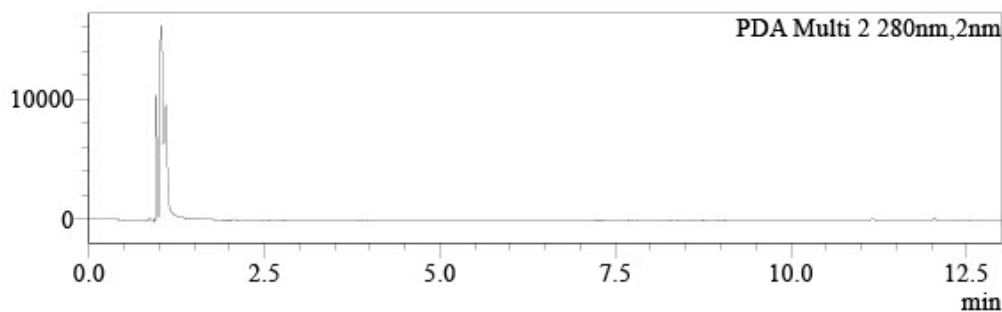

Peak Table

PDA Ch1 225nm

Name	Ret. Time	Area	Conc.	Unit	Area%
GHK-Cu	1.017	5595520	52.895	mg	79.795
	1.326	2284	0.000		0.033
	1.551	1219	0.000		0.017
	1.820	150	0.000		0.002
	2.076	1991	0.000		0.028
	2.568	255	0.000		0.004
	2.773	1352	0.000		0.019
	8.431	217	0.000		0.003
	10.133	635	0.000		0.009
	10.468	214	0.000		0.003
	10.767	168	0.000		0.002
	10.884	1848	0.000		0.026
BPC-157	11.159	851378	13.359	mg	12.141
	11.414	2485	0.000		0.035
	11.595	4914	0.000		0.070
	11.824	581	0.000		0.008
	11.936	247	0.000		0.004

Name	Ret. Time	Area	Conc.	Unit	Area%
TB-500	12.034	544770	11.696	mg	7.769
	12.201	345	0.000		0.005
	12.534	1251	0.000		0.018
	12.909	559	0.000		0.008
		7012383			100.000

Peak Table

PDA Ch2 280nm

Name	Ret. Time	Area	Conc.	Unit

Analysis Report




LIQUILABS

Analytical services


Analysis of identity of active ingredient by UHPLC with mass spectrometric detection

MS Chromatogram

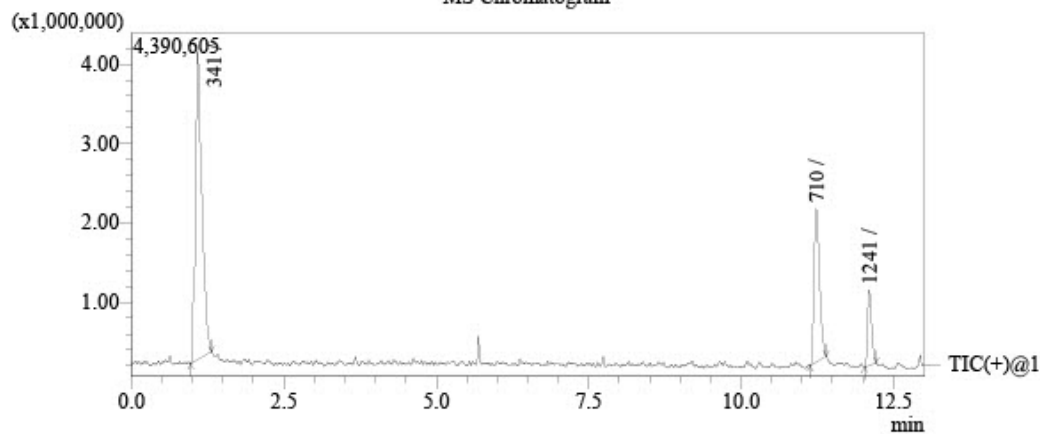




Library Search



Detected compound
Ret. Time: 0.974 Base Peak: 341
Scan Range: 280.00 - 2000.00 Polarity: Positive

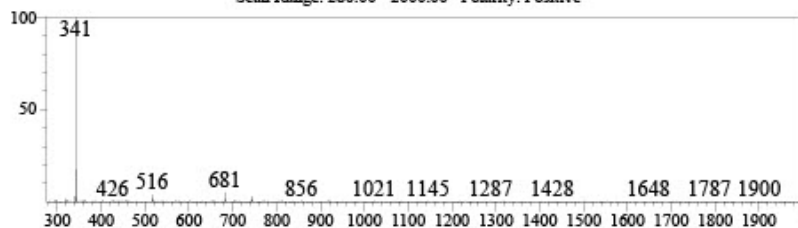

Library: Peptide screening lib
Formula: C₁₄H₂₄CuN₆O₄ CAS: 49557-75-7 Mol.Weight: 401 Ret.Index: 0
Compound Name: GHK-Cu

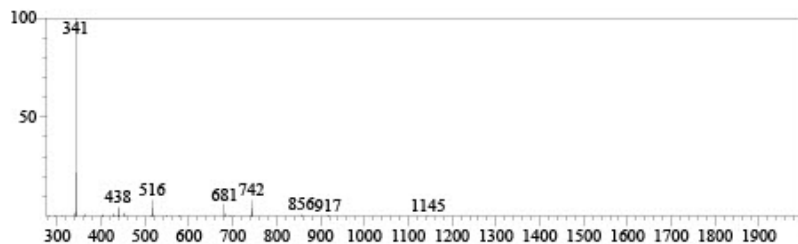

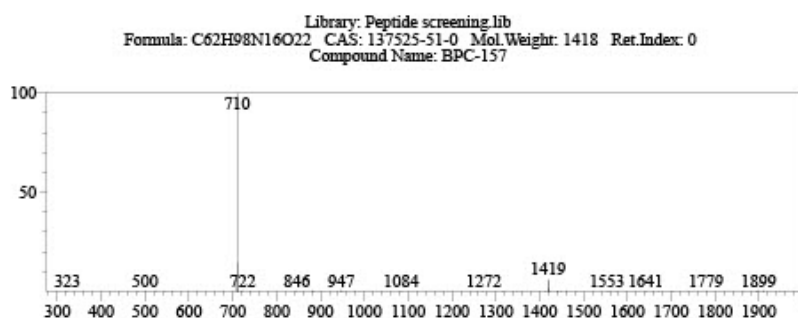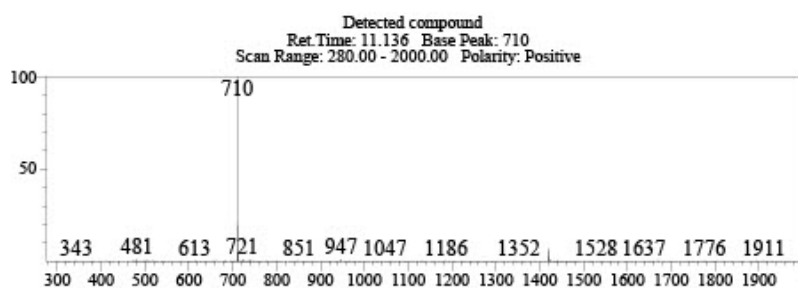

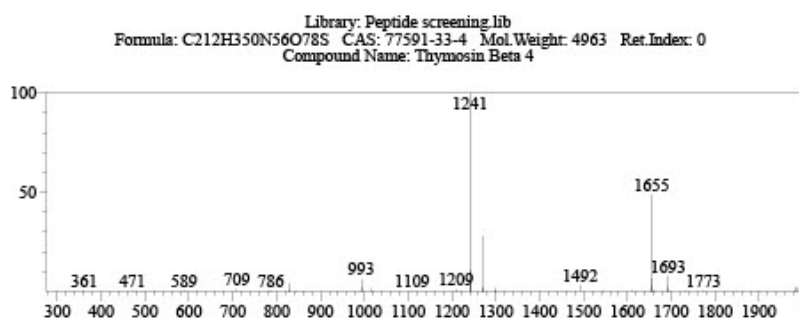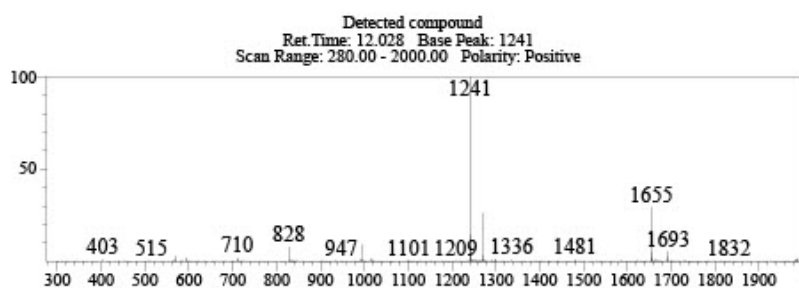

 LIQUILABS Analytical services	Method Specification	
Determination of bioburden of lyophilized samples		
Document number MIC_001_2025	Superseded document -	Number of pages 2

1. Instrumentation and chemicals

1.1. Instruments used

- Sterile Syringe 2mL Luer
- Sterile needles
- Ready made PCA Plate ROTI Aquatest
- Ready made Sab4 Plate ROTI Aquatest

1.2. Chemicals

Sterile physiological solution (0.9% NaCl)

2. Sample preparation and inoculation

2.1 Sample preparation

1. Fresh sterile needle and syringe was used for measuring exactly 2 mL of sterile physiological solution.
2. Needle was changed and by new needle rubber top of peptide container was penetrated and 2 mL of sterile physiological solution was dispensed.
3. Content of container was completely dissolved and left for 5 minutes to settle potentially created bubbles.
4. This procedure is repeated for two vials.

2.2 Total Aerobic microbial count inoculation and cultivation

1. By sterile needle 1 mL of solution was filled into the sterile syringe.
2. Needle was placed above the flame for few seconds to sterilize.
3. Consequently 1 mL of solution was poured into the ready to use sterile petri dish filled with PCA agar and petri dish was closed.
4. Process was repeated for two petri dishes.
5. With sterile needle, 1 mL of sterile physiological solution was filled into the sterile needle and was inoculated onto one sterile petri dish filled with PCA agar as negative control sample.
6. Samples and negative control sample were placed in incubator at temperature 37°C for 120h.

2.3 Total Yeast and Mold count inoculation and cultivation

1. By sterile needle 1 mL of solution was filled into the sterile syringe.
2. Needle was placed above the flame for few seconds to sterilize.
3. Consequently 1 mL of solution was poured into the ready to use sterile petri dish filled with Sab4 agar and petri dish was closed.
4. Proces was repeated for two petri dishes.
5. With sterile needle, 1 mL of sterile physiological solution was filled into the sterile needle and was inoculated onto one sterile petri dish filled with Sab4 agar as negative control sample.
6. Samples and negative control sample were placed in incubator at temperature 25°C for 72h.

3. Evaluation of results

After incubation time, colonies are counted as cfu (colonies forming units) and result per 1 g of sample is determined as:

$$CFU_{avg} = \frac{\sum CFU_n}{n}$$

CFU_{avg} = average CFU counted form n inoculations

CFU_n = CFU counted per inoculation

n = number of inoculations

$$CFU \text{ per gram} = \frac{CFU_{avg}}{m_s} * DF$$

CFU_{avg} = Average CFU counted from n inoculations

m_s = mass of sample (mg)

DF = Dilution factor

If negative control sample is evaluated as positive, process have to be repeated due to possible contamination in the process of inoculation or incubation.

Responsibles



Mr. Ján Galbavý
CEO

Analysis results relate only to the samples tested.

This document shall not be reproduced except in full, without the written approval of Liquilabs s.r.o.